**MULTI SCIENCES**

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**Mouse Th17 Staining Kit**

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This package insert must be read entirely before using this product. For proper performance, use the insert provided with each individual product received.

**Catalog Number**

KTH217 - 25

KTH217 - 100

Optimization for mouse Th17 staining in anticoagulated blood and splenocytes.

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**For research use only. Not for use in diagnostic procedures.**

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## INTRODUCTION

The kit is designed to detect frequencies of Th17 cells in T lymphocytes (or CD4<sup>+</sup> helper T cells) in mouse anticoagulated blood or splenocytes. Th1, Th2 and Th17 cells are subsets differentiated from helper T cells (strictly speaking, that is Th0 cells) under physiological and pathological conditions. There are extremely few Th1, Th2 or Th17 cells in peripheral blood in a resting state (unstimulated state, e.g. normal physiological state in mouse) due to the weak differentiation ability of Th0 cells, thus, it's hardly to detect IFN- $\gamma$ , IL-4 and IL-17A. Th0 cells could differentiate to Th1, Th2 or Th17 cells dependent on different cytokines when stimulated by external factors (e.g. stimuli, pathogens), and more IFN- $\gamma$ , IL-4 and IL-17A could be detected. In experiment, it is to measure the ability of Th cells in response to stimuli.

Generally, Phorbol 12-Myristate 13-Acetate (PMA) and Ionomycin are used as stimuli. PMA, a PKC activator, could activate phosphorylation of many protein kinases in downstream, cascade reaction induces protein expression, then results in activation of T cells. PKC could be activated in the combination of DAG and Ca<sup>2+</sup> in cells, with the involvement of Ionomycin (a calcium transporter, transport Ca<sup>2+</sup> from organelles to cytoplasm), PKC could be further activated in T cells. So PMA and Ionomycin could activate T cells synergistically.

Activated T cells would secrete cytokines to extracellular environment, which could not be detected by flowcytometry, thus cytokines should be blocked in the cell. Cytokines are synthesized in the Golgi and some may traffic through the endoplasmic reticulum to be released in soluble form in the extracellular environment. Blocking secretion by disrupting Golgi to cut off the way of cytokine transport. Generally, Brefeldin A and/or Monensin are used as blocker.


Sometimes, more cytokines could be detected in sample after separation by lymphocyte separation medium in our experiment. The best protocol should be optimized by the end-user.

## LIMITATIONS OF THE PRODUCT

The product is intended for flowcytometry applications. FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES.

The product should not be used beyond the expiration date on the label.

Do not mix reagents with those from other lots or sources.

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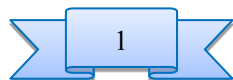
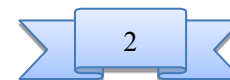
**Note:** PMA/Ionomycin Mixture (250 $\times$ ) and BFA/Monensin Mixture (250 $\times$ ) can be stored at - 20 $^{\circ}$ C, FIX & PERM Medium A and Medium B can be stored at room temperature, other reagents can be stored at 2 - 8 $^{\circ}$ C. All reagents can be stable for at least 1 year when stored at recommended condition.

### OTHER SUPPLIES REQUIRED

1. **Cell culture medium** (e.g. RPMI 1640, DMEM).
2. **Fetal Bovine Serum** (Maybe needed)
3. **Lymphocyte separation medium** (Cat No: MLSM1092, MultiSciences) (Maybe needed)
4. **Compensation Beads** (Maybe needed)
5. **Paraformaldehyde, 4%** (Cat No: F0001, MultiSciences) (Maybe needed)
6. **12  $\times$  75 mm round bottom test tubes.**
7. **Vortexer.**
8. **CO<sub>2</sub> incubator.**
9. **Swing-out horizon centrifuge** (with rotor for 15 ml tubes).

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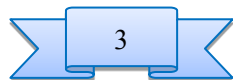## PROTOCOL

### Sample collection and treatment

Collect appropriate amount of anticoagulant blood by using heparin as anticoagulant or other anticoagulants (e.g. EDTA, sodium citrate), store at room temperature or 2 - 8°C, and detect on the same day.

For spleen, collect and prepare as single cell suspension by using appropriate method, and remove tissue mass, detect on the same day.

### Sample preparation

1a. *For anticoagulant blood with heparin as anticoagulant*, pipet 125  $\mu$ l anticoagulant blood into 12  $\times$  75 mm round bottom test tubes, add 125  $\mu$ l medium without FBS and 1  $\mu$ l PMA/Ionomycin Mixture (250 $\times$ ) and 1  $\mu$ l BFA/Monensin Mixture (250 $\times$ ). Mix 125  $\mu$ l anticoagulant blood and 125  $\mu$ l medium, served as control. Mix well, incubate in CO<sub>2</sub> incubator at 37°C for 4 - 6 hours, vortex every 1 - 2 hours during incubation.


**Note:** Anticoagulant blood with heparin as anticoagulant could also be treated by method 1b, more cytokines could be detected in our tests. The best protocol should be optimized by the end-user.

1b. *For anticoagulant blood with other anticoagulant (e.g. EDTA, sodium citrate)*, isolate peripheral blood mononuclear cells (PBMCs) by using Lymphocyte separation medium (Cat No: MLSM1092, MultiSciences). Resuspend pellet at 1  $\times$  10<sup>7</sup>/ml in medium with 10 % FBS. Pipet 250  $\mu$ l PBMCs into 12  $\times$  75 mm round bottom test tubes, add 1  $\mu$ l PMA/Ionomycin Mixture (250 $\times$ ) and 1  $\mu$ l BFA/Monensin Mixture (250 $\times$ ). With PBMCs only as control. Mix well, incubate in CO<sub>2</sub> incubator at 37°C for 4 - 6 hours, vortex every 1 - 2 hours during incubation.

1c. *For spleen*, prepare single cell suspension by using appropriate method and remove tissue mass. (Optional) Isolate spleen mononuclear cells by using Lymphocyte separation medium (Cat No: MLSM1092, MultiSciences).

Resuspend pellet at 1  $\times$  10<sup>7</sup>/ml in medium with 10 % FBS. Pipet 250  $\mu$ l mononuclear cells into 12  $\times$  75 mm round bottom test tubes, add 1  $\mu$ l PMA/Ionomycin Mixture (250 $\times$ ) and 1  $\mu$ l BFA/Monensin Mixture (250 $\times$ ). With mononuclear cells only as control. Mix well, incubate in CO<sub>2</sub> incubator at 37°C for 4 - 6 hours, vortex every 1 - 2 hours during incubation.

**Note:** Sometimes, more cytokines could be detected in mononuclear cells after isolation by Lymphocyte separation medium in our tests. The best protocol should be optimized by the end-user. PMA/Ionomycin Mixture (250 $\times$ ) and BFA/Monensin Mixture (250 $\times$ ) are extremely volatile, tighten immediately the lid after use.

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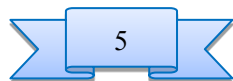
- Pipet 100  $\mu$ l sample or control into new test tubes, add 5  $\mu$ l Anti-Mouse CD3 $\epsilon$ , FITC and 5  $\mu$ l Anti-Mouse CD4, PerCP-Cy5.5 to each tube. Vortex to mix well, and incubate at room temperature for 15 minutes, protect from light.
- Add 100  $\mu$ l FIX & PERM Medium A to each tube. Vortex to mix well, and incubate at room temperature for 15 minutes, protect from light.
- Add 3 ml pre-cool 1 $\times$  Flow Cytometry Staining Buffer to each tube. Centrifuge at 300  $\times$  g for 5 minutes, discard supernant.  
**Note:** Discard supernant as much as possible, with little residual liquid left.
- Add 100  $\mu$ l FIX & PERM Medium B and 5  $\mu$ l Anti-Mouse IL-17A, PE to each test tube. Vortex to mix well, and incubate at room temperature for 15 minutes, protect from light.
- Add 2 ml 1 $\times$  Flow Cytometry Staining Buffer to each tube. Centrifuge at 300  $\times$  g for 5 minutes, discard supernant.
- Resuspend pellet by adding 500  $\mu$ l 1 $\times$  Flow Cytometry Staining Buffer to each tube, ready for determination. Or add 500  $\mu$ l 1 - 4 % paraformaldehyde to resuspend pellet, then store at 2 - 8 $^{\circ}$ C, protect from light, determine within 24 hours.

#### Detection by flowcytometry

- Voltage and compensation adjustment  
 It is recommended to adjust voltage and compensation of flow cytometer by using compensation beads. It is important to adjust voltage of FSC/SSC and gate target cells because compensation beads are often smaller than cells.
- Correct gating to determine the frequencies of Th17 cells in CD3 $^{+}$ CD4 $^{+}$  T cells.  
**Note:** Acquire at least 20,000 to 30,000 CD3 $^{+}$ CD4 $^{+}$  T cells. Depending on the species and individual differences, frequencies of cytokine producing cells derived from activation of mouse lymphocytes can vary widely for a particular cytokine. In order to make statistically significant frequency measurements, sufficiently large sample sizes should be acquired during flow cytometric analysis.

#### Excitation wavelength and emission wavelength of fluorophore

Fluorophore	Ex (nm)	Em (nm)
FITC	495	519
R-Phycoerythrin (PE)	480;565	578
PerCP-Cy5.5	488	682

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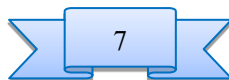
Examples of results/结果示例

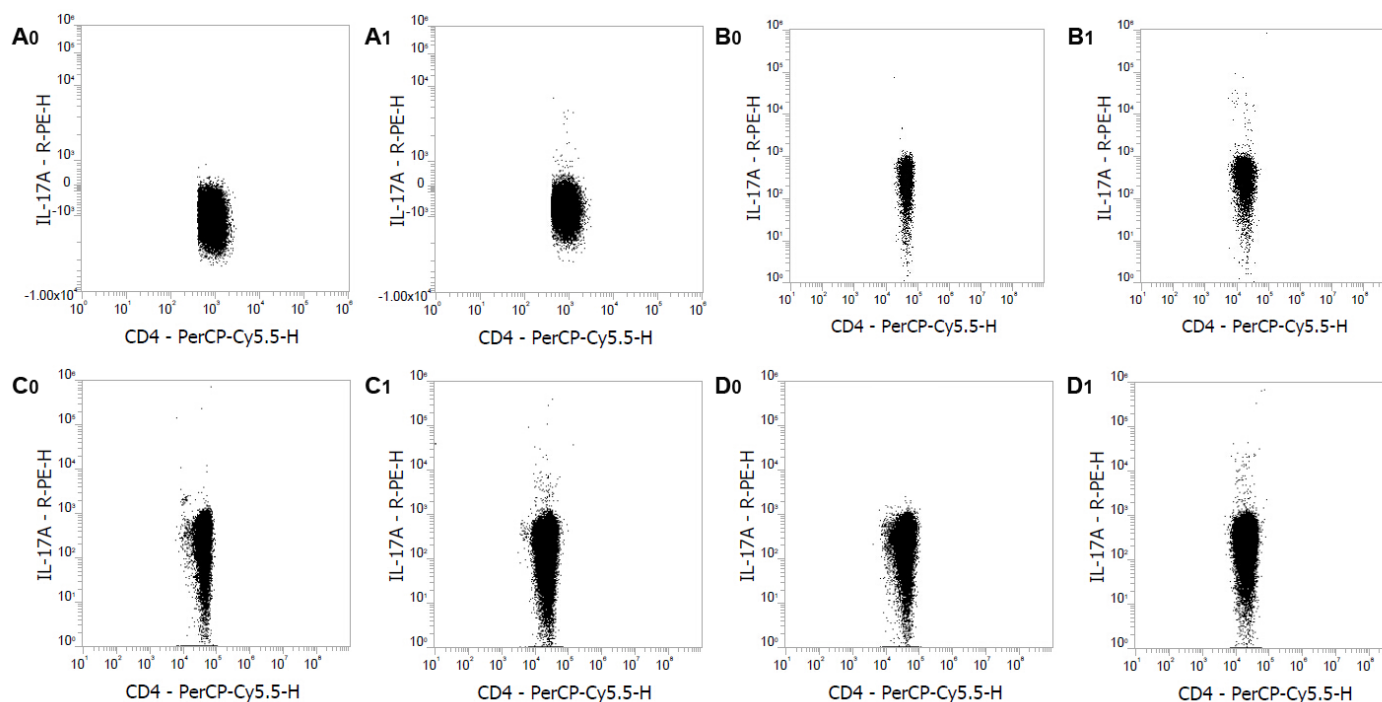**Flow cytometric analysis of Mouse Th17 Staining Kit.** The staining pattern of IL-17A on resting heparin anticoagulated blood (A<sub>0</sub>), PBMCs from EDTA anticoagulated blood (B<sub>0</sub>), splenocytes (C<sub>0</sub>) and spleen mononuclear cells (D<sub>0</sub>) of normal ICR mouse. PMA/Ionomycin stimulated heparin anticoagulated blood (A<sub>1</sub>), PBMCs from EDTA anticoagulated blood (B<sub>1</sub>), splenocytes (C<sub>1</sub>) and spleen mononuclear cells (D<sub>1</sub>) of normal ICR mouse. Dot plot analysis are derived from gated CD3<sup>+</sup>/CD4<sup>+</sup> cells population. Flow cytometry was performed on a Thermo Fisher Attune NxT.

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